


SlimWhite

Bright and clear

- **Uniform:**
Illumination is designed by light guide plate with diffuser
- **Easy to observe:**
Light source is under the gel, avoid bubble and reflection
- **Cost Effective:**
Long lifetime (~50,000 hrs), low power consumption



SlimWhite white light panel, the light source which can be used independently, easy to view colorimetric gels, coomassie blue or silver stained media. No need to use UV for white light conversion.

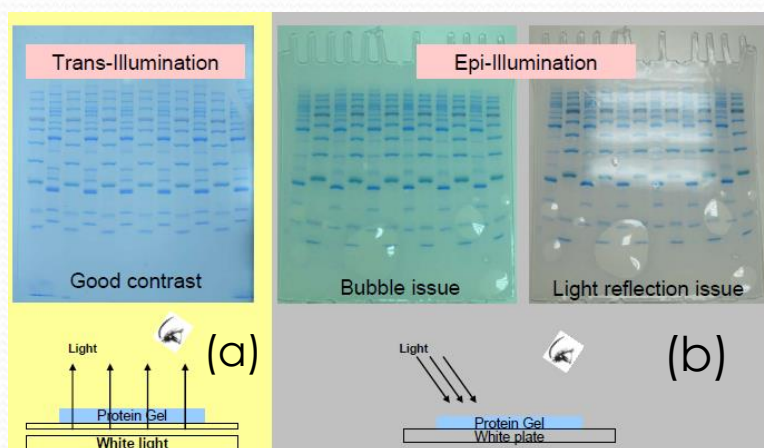

- (a) Light source is under the gel, avoid bubble and reflection.
(b) Ordinary light source is upon the gel with bubble and reflection affect.

Order Information and Specification:

Products	SlimWhite
Cat No.	SW-01
Outer dimensions (mm)	210D x 305W x 15H
Illumination area (mm)	170D x 275W
Light source	LED s
Power	12V/DC, 0.5A
Weight (Kg)	1.3 kg
Wavelength (nm)	White (blue chip + yellow fluorescent powders)
Battery box	Optional