

SlimWhite

ILLUMINATOR



A whole new design of SlimWhite white light panel, the light source which can be used independently, easy to view colorimetric gels, coomassie blue or silver stained media. No need to use UV for white light conversion. High brightness and low power consumption.

● **Uniform:**

Illumination is designed by light guide plate with diffuser


● **Easy to observe:**

Light source is under the gel, avoid bubble and reflection issue

● **Cost Effective:**

Long lifetime (~50,000 hrs), low power consumption


Trans-illumination

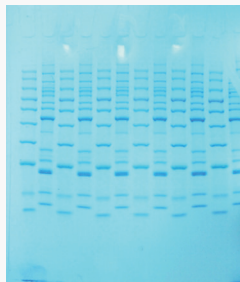

Good contrast


(O)

Epi-illumination

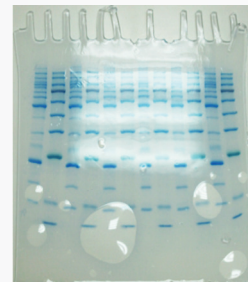


Light reflection issue

(X)

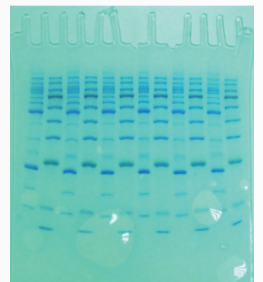

Bubble issue

(X)

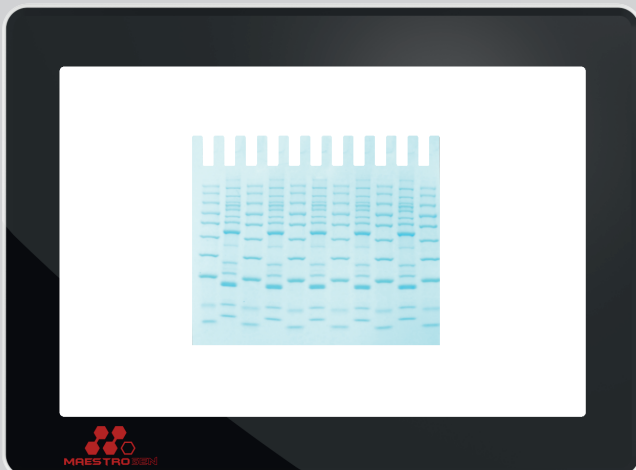

Specification

Cat no.: SW-02

Dimensions (mm): 292 (D)* 217 (W)* 12.8 (H)

View Size (mm): 240 (D)* 160 (W)

Light Source: White light LEDs

Power: DC12V

Weight (Kg): 1.2